

Quantification of Chlorophyll as a Proxy for Biofilm Formation in the Cyanobacterium *Synechococcus elongatus*

Eleonora Sendersky¹, Ryan Simkovsky², Susan S. Golden² and Rakefet Schwarz^{1,*}

¹The Mina and Everard Goodman Faculty of Life Sciences, Bar-Ilan University, Ramat-Gan, Israel;

²Division of Biological Sciences, University of California, San Diego, La Jolla, CA, USA

*For correspondence: Rakefet.Schwarz@biu.ac.il

[Abstract] A self-suppression mechanism of biofilm development in the cyanobacterium *Synechococcus elongatus* PCC 7942 was recently reported. These studies required quantification of biofilms formed by mutants impaired in the biofilm-inhibitory process. Here we describe in detail the use of chlorophyll measurements as a proxy for biomass accumulation in sessile and planktonic cells of biofilm-forming strains. These measurements allow quantification of the total biomass as estimated by chlorophyll level and representation of the extent of biofilm formation by depicting the relative fraction of chlorophyll in planktonic cells.

Keywords: Biofilm, Cyanobacteria, *Synechococcus elongatus*, Chlorophyll measurement, Sessile, Planktonic

[Background] Several recently published studies indicate an emerging interest in the mechanisms that underlie cell-aggregation and biofilm development in cyanobacteria (Fisher *et al.*, 2013; Jittawuttipoka *et al.*, 2013; Schatz *et al.*, 2013; Enomoto *et al.*, 2014; Schwarzkopf *et al.*, 2014; Enomoto *et al.*, 2015; Oliveira *et al.*, 2015; Agostoni *et al.*, 2016; Parnasa *et al.*, 2016). We recently reported a self-biofilm-inhibitory mechanism that dictates planktonic growth of the model unicellular cyanobacterium *Synechococcus elongatus* PCC 7942 (Schatz *et al.*, 2013; Nagar and Schwarz, 2015). Abrogation of the biofilm-inhibitory process by inactivation of particular genes results in robust biofilm development in this otherwise planktonic strain (Schatz *et al.*, 2013; Nagar and Schwarz, 2015). These studies required quantification of the extent of biofilm development in various strains and under different conditions. Crystal violet is commonly used for quantification of biofilms in heterotrophic bacteria (O'Toole and Kolter, 1998). This staining procedure, however, quantifies only the sessile fraction of cells. Here we provide a detailed protocol for culture growth and quantification of cyanobacterial biofilms using chlorophyll measurement as a proxy for biomass accumulation in sessile as well as in planktonic cells. These measurements allow estimation of the total biomass accumulated and representation of the relative fraction of chlorophyll in sessile or in planktonic cells.

Materials and Reagents

1. Custom-made Pyrex glass tubes for bacterial liquid cultures (200 x 32 mm, made from Pyrex tubing, Corning, catalog number: 8510-32-D)

Note: Can order from <http://www.degroot.co.il/>.

2. 'Sponge plug' (Plastic foam stoppers, 27 x 34 mm) (Jaece Industries, catalog number: L800-C)
3. Pasteur pipettes (230 mm) (Romical, catalog number: 94-08401002)
4. 0.45 µm syringe filter (Sartorius, catalog number: 16555)
5. Filter Stericup-GP 250 ml Express Plus PES (0.22 µm) (EMD Millipore, catalog number: SCGPU05RE)
6. Silicone tubing (6 x 9 mm, 4 x 6 mm) (Degania Silicone, catalog numbers: 2110600234, 2110400434, respectively)
7. Sterile pipettes 1 and 25 ml
 - 1 ml pipettes (Corning, Costar[®], catalog number: 4011)
 - 25 ml pipettes (Corning, Costar[®], catalog number: 4489)
8. Sterilized pipette tips (Corning, Axygen[®], catalog numbers: T-200-C, T-1000-C)
9. Eppendorf tubes (1.5 ml) (Corning, Axygen[®], catalog number: MCT-175-C)
10. *Synechococcus elongatus* PCC7942
11. Sodium nitrate (NaNO₃) (Sigma-Aldrich, catalog number: S8170)
12. Magnesium sulfate heptahydrate (MgSO₄·7H₂O) (Merck, catalog number: K26364082)

Note: This product has been discontinued. Alternatively, Merck, catalog number: 105886 can be used.
13. Calcium chloride dihydrate (CaCl₂·2H₂O) (ICN, catalog number: 10035-04-8)

Note: This product has been discontinued. Alternatively, Bio-Lab, catalog number: 034205 can be used.
14. Potassium phosphate dibasic (K₂HPO₄) (Honeywell International, Riedel-de-Haen, catalog number: 04248)
15. Ethylenediaminetetraacetic acid (Na₂Mg-EDTA) (Bio-Lab, catalog number: 05142359)

Note: This product has been discontinued. Alternatively, Biosolve, catalog number: 051423 can be used.
16. Ferric ammonium citrate (C₆H₁₁FeNO₇) (MP Biomedicals, catalog number: 02158040)
17. Citric acid (C₆H₈O₇) (Frutarom, catalog number: 2355511000)

Note: This product has been discontinued. Alternatively, Biosolve, catalog number: 030205 can be used.
18. Boric acid (H₃BO₃) (Bio-Lab, catalog number: 02010591)

Note: This product has been discontinued. Alternatively, Biosolve, catalog number: 020105 can be used.
19. Manganese chloride hexahydrate (MnCl₂·6H₂O) (Duchefa Biochemie, catalog number: M0533)
20. Zinc sulfate heptahydrate (ZnSO₄·7H₂O) (CARLO ERBA Reagents, catalog number: 494907)
21. Sodium molybdate dihydrate (Na₂MoO₄·2H₂O) (MP Biomedicals, catalog number: 194863)
22. Copper sulphate pentahydrate (CuSO₄·5H₂O) (Honeywell International, Riedel-de-Haen, catalog number: 12849)

23. Cobalt(II) nitrate hexahydrate ($\text{Co}(\text{NO}_3)_2 \cdot 6\text{H}_2\text{O}$) (Honeywell International, Riedel-de-Haen, catalog number: 12922)
Note: This product has been discontinued. Alternatively, Sigma-Aldrich, catalog number: 239267 can be used.
24. N-2-hydroxyethylpiperazine-N'-2-ethanesulfonic acid (HEPES) (Sigma-Aldrich, catalog number: H3375)
25. Sodium hydroxide (NaOH) (Frutarom, catalog number: 5553510)
26. Acetone ($(\text{CH}_3)_2\text{CO}$), A.C.S. reagent (Avantor Performance Materials, J.T. Baker®, catalog number: 9006-03)
27. Sodium bicarbonate (NaHCO_3) (DAEJUNG CHEMICAL & METALS, catalog number: 7566-4100)
28. Sodium thiosulfate ($\text{Na}_2\text{S}_2\text{O}_3$) (Sigma-Aldrich, catalog number: 217247)
29. Bacto agar (BD, Bacto™, catalog number: 214010)
30. LB agar (Lennox) (BD, Difco™, catalog number: 240110)
31. Antibiotics (added as appropriate according to the resistance of the particular strain)
 - a. Spectinomycin dihydrochloride pentahydrate (Duchefa Biochemie, catalog number: S0188)
 - b. Kanamycine sulphate monohydrate (Duchefa Biochemie, catalog number: K0126)
 - c. Gentamycin sulphate (Duchefa Biochemie, catalog number: G0124)
 - d. Chloramphenicol (Duchefa Biochemie, catalog number: C0113)
32. Liquid BG11-medium (see Recipes)
 - a. Stock I (100x concentrated)
 - b. Stock II (100x concentrated)
 - c. Stock III (100x concentrated)
 - d. Stock V (1,000x concentrated)
33. Solid BG11-medium (see Recipes)

Equipment

1. Laminar flow hood
2. Fume hood
3. Autoclave
4. Spectrophotometer (Agilent Technologies, model: Cary 100, catalog number: 10069000) and respective cuvettes (Cell type P.L = 10 mm/EA) (Starna Cells, catalog number: 9-SOG-10)
5. Benchtop centrifuge (MiniSpin, max. centrifugal force: 12,100 x g) (Eppendorf, model: MiniSpin®, catalog number: 5452000018)
6. Refrigerator (4-7 °C)
7. Growth rooms (30 ± 2 °C and 24 ± 2 °C for liquid cultures and cultures on solid medium, respectively)

Note: Growing cultures on solid medium at 24 °C rather than 30 °C allows maintaining the cultures for longer periods.

8. Light source

For liquid cultures—incandescent light producing a flux of 20-30 $\mu\text{mol photons m}^{-2} \text{sec}^{-1}$; Cultures on solid-BG11 are illuminated by fluorescent light (5 $\mu\text{mol photons m}^{-2} \text{sec}^{-1}$)

9. Pipet-aid

10. Single-channel pipettes 200 and 1,000 μl

11. Bubbling system

- a. Air compressor (Oil free air compressor) (Assouline Compressors, model: vs 204 50)
- b. Cylinder carbon dioxide (CO_2 , compr. 99.5%) (Maxima, catalog number: GCDU2.527)
- c. Flow meters (Tuttnauer company, 0-0.1 L/min, 1-10 L/min)
- d. High-flow air pressure regulator (0-2 PSI) (Marsh Bellofram, model: Type 70, catalog number: 960-129-000)

Note: To bubble cultures with 3% CO_2 in air the following setup is used (Figure 1): CO_2 gas is mixed with compressed air using flow meters to yield 3% CO_2 in air. This mixture is humidified by bubbling into 1.5 L double distilled water (DDW) in a 2 L bottle to prevent water evaporation during bubbling of the cultures. The outlet from this bottle is passed through an empty 2 L bottle (serving to trap residual liquid) and through an air pressure regulator (200 mbar is appropriate for bubbling of ~300 culture tubes). Bubbling of CO_2 enriched air from the latter into multiple cultures is obtained using home-made manifolds constructed from silicone tubing and appropriate T-bar connectors. The gas in each manifold channel is passed through a 0.45 μm filter (Figure 1). A short video is provided to demonstrate the flow rate in individual cultures tubes (Video 1).

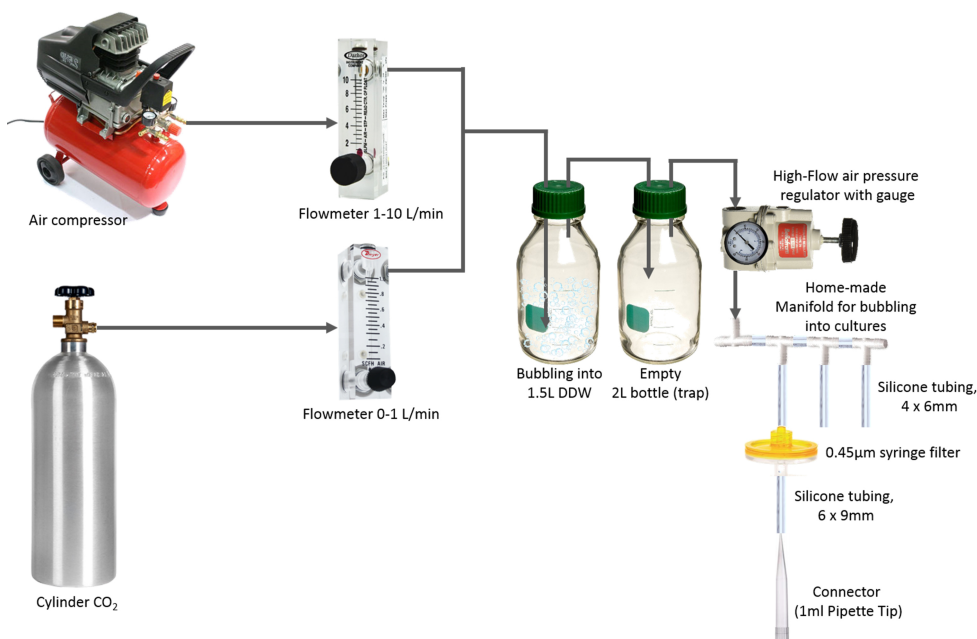


Figure 1. A scheme describing the setup used to bubble cultures with CO_2 -enriched air

Video 1. Cyanobacterial cultures under bubbling



Procedure

1. Growing liquid cultures of *Synechococcus elongatus* PCC 7942 and biofilm-forming mutants thereof
 - a. Round-bottom Pyrex tubes (20 cm in length, 3 cm in diameter) are used as growth vessels.
 - b. A cotton-plugged Pasteur pipette (23 cm long) is inserted through a sponge-plug such that the tip of the Pasteur pipette is placed 1-2 cm above the bottom of the tube (Figure 2A, 'growth set').
 - c. Growth sets are autoclaved in liquid cycle.
 - d. Sterile liquid BG11-medium (25-50 ml) (see Recipes) is transferred into the Pyrex tube. The culture volume is adjusted to the biomass required for the designated experiment.
 - e. Cells are scraped from solid-BG11 medium (see Recipes) with a 1 ml sterile plastic pipette and the biomass is shaken in the liquid medium (see Note 4).

Important: To reproducibly observe biofilm formation, it is crucial to add the ferric ammonium citrate and citric acid components from a freshly made stock (see note in Recipe 1f). Autoclaved BG11 should be used within 4 days to inoculate cultures.

- f. Connect the growth set to the bubbling system (see description above) via 1 ml pipette tip (Figure 2A).
- g. The use of manifolds allows simultaneous bubbling of multiple cultures (Figure 2B, see Note under Equipment).

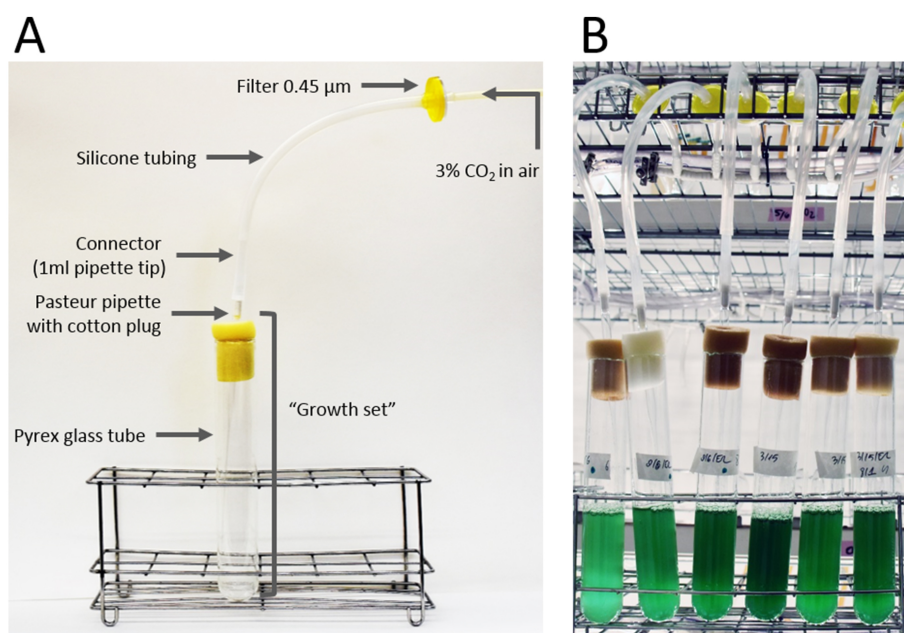


Figure 2. Growing cyanobacterial cultures. A. Growth set connected to a single manifold channel; B. Cultures in the growth room.

2. Quantification of biofilm

- a. Cells are cultured under continuous bubbling of 3% CO₂ in air at 30 °C under continuous illumination with incandescent light (20-30 µmol photons m⁻² sec⁻¹).
- b. Starting cultures at the exponential phase are diluted (usually following 3 days of growth) to an optical density at 750 nm of 0.5.
- c. Experiments are initiated by an additional dilution of the cultures to an optical density at 750 nm of 0.5 (see Note 5).
- d. Biofilm development under this setup typically begins after 2-3 days.
- e. To verify that cultures are axenic 'spot' a 20 µl aliquot on solid LB agar.
- f. Percentage of chlorophyll in suspended cells serves to quantify biofilm formation as follows (Figure 3):
 - i. Extraction of chlorophyll in planktonic cells is performed in 1.5 ml Eppendorf tubes. Remove an aliquot from the upper part of the suspended culture (see Note 6). In cases where the suspended fraction appears especially dense (planktonic strains or poor biofilm formers), a 0.2 ml sample is used for extraction. When robust biofilms are formed, 15 ml of planktonic cells are removed and concentrated 3 to 6-fold by centrifugation (5,000 x g, 10 min) prior to the chlorophyll extraction in 80% acetone (final concentration).
 - ii. Extraction of chlorophyll in the biofilm is performed in the growth tube. To separate the biofilm fraction, slowly remove the suspended cells. Biofilms tend to collapse during removal of the last few ml of liquid. Therefore, once clumps of biofilm detach and fall into the liquid, stop removing the culture. Measure the volume remaining in the tubes

using a pipette (in our hands 2-5 ml are left) as well as the volume removed from the culture. The sum of these volumes (V_{total}) serves for calculation of the amount of chlorophyll in the suspended cells (see Data analysis). Extraction of chlorophyll is performed in 80% acetone at a final volume equal to V_{total} .

- iii. Extraction of chlorophyll is carried out overnight in the refrigerator and chlorophyll is quantified based on absorbance as previously described (Arnon, 1949; Ritchie, 2006). Chlorophyll may be extracted with methanol instead of acetone. A step-by-step protocol for measurement of chlorophyll a concentration in cyanobacteria is described by Zavřel *et al.* (2015).

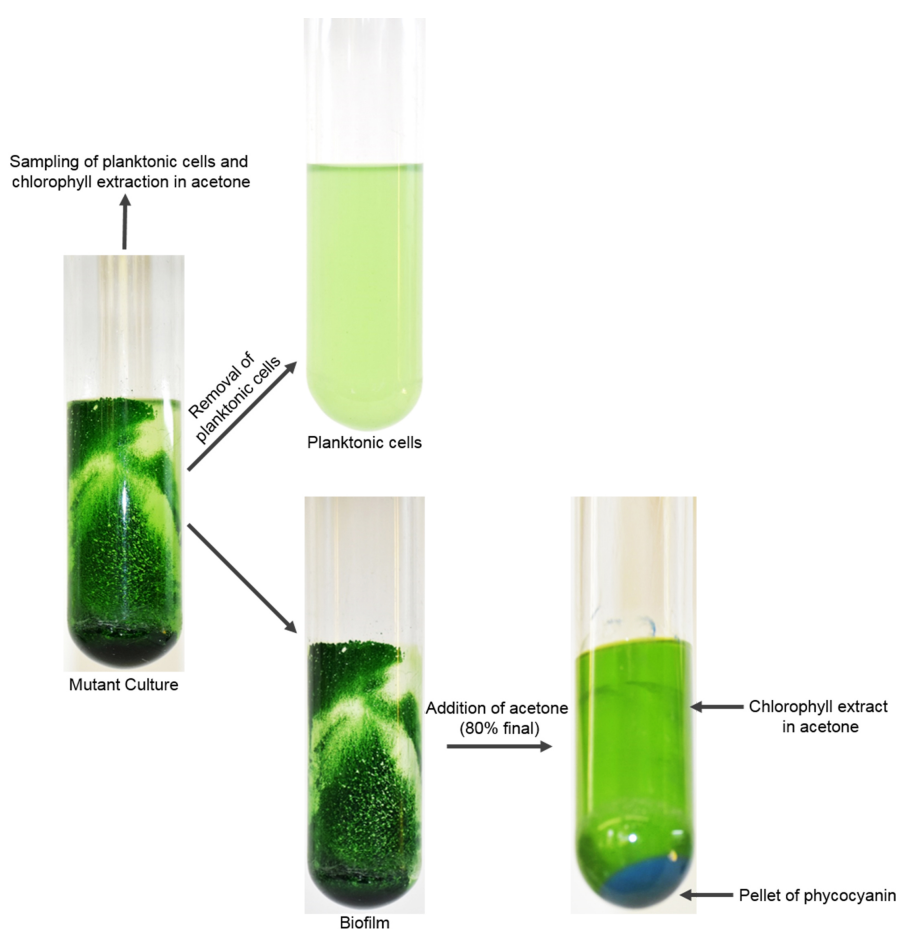


Figure 3. Quantification of biofilm based on chlorophyll measurement

Data analysis

The percentage of chlorophyll in suspension is calculated from the amount of chlorophyll in the planktonic and biofilm fractions. These values are calculated based on measured chlorophyll absorbances, sample volumes, and dilution factors, as described below. An example calculation is provided in Table 1.

Table 1. An example of calculations of chlorophyll (Chl) in planktonic and biofilm cells

Strain	Fraction	Volume (ml)			Extraction in Acetone			Dilution of extract ^a		Dilution factor ^b	OD ₆₆₃	OD ₇₅₀	Chl concentration (µg/ml) ^c	Total Chl (µg) ^d	Chl in Suspension (%)
		Suspension	Detached biofilm	Total Volume	Planktonic Aliquot (ml)	Added DDW (ml)	100% acetone added (ml)	Aliquot from Chl Extract (ml)	80% acetone added (ml)						
Wild type [#]	Planktonic	50		50	0.05	0.15	0.8			20	0.1257	0.0013	31.60	1579.88	100.00
Biofilm forming mutant - without Planktonic Correction	Planktonic	48		50	0.2	0	0.8			5	0.0412	0.0016	2.51	125.73	8.49
	Biofilm		2			8	40	0.1	0.9	10	0.2147	0.0013	27.10	1355.09	
Biofilm forming mutant - Planktonic Corrected	Planktonic	48		50	0.2	0	0.8			5	0.0412	0.0016	2.51	125.73	8.52
	Biofilm		2			8	40	0.1	0.9	10	0.2147	0.0013	27.10	1350.06	

[#]Grows planktonically.

^aDilution of Chl extract is performed to reach the linear range of the assay.

^bDilution factor is calculated from the dilution of the aliquot during acetone extraction and the dilution of the Chl extract for spectrophotometry.

^cChl concentration = (OD₆₆₃ - OD₇₅₀) x 12.7 x (Dilution Factor); see Arnon, 1949.

^dFor planktonic samples, (Total Chl) = (Chl concentration) x (Total Volume). For biofilm samples from robust strains, the contribution of Chl from planktonic cells is negligible, and thus (Total Chl) = (Chl concentration) x (Volume of Chl extract), where 'Volume of Chl extract' should be equivalent to 'Total Volume'. If correcting for the contribution of planktonic Chl to the biofilm fraction, then (Total Chl) = (Chl concentration in biofilm) x (Volume of Chl extract) - (Chl concentration in planktonic) x (Volume of Detached biofilm).

1. Calculate the amount of chlorophyll in suspended cells based on measured chlorophyll concentration and the volume of the suspended culture (V_{total}).

$$\text{Total Chl Planktonic } (\mu\text{g}) = (\text{OD}_{663} - \text{OD}_{750}) \times 12.7 \times (\text{Dilution Factor}) \times V_{total}$$

2. Calculate the amount of chlorophyll in the biofilm based on measured chlorophyll concentration and the volume of the acetone extraction, which should be V_{total} .
 - a. If robust biofilm formation is observed, then the contribution from the remaining planktonic cells in the biofilm fraction is negligible and the total chlorophyll is calculated as follows:

$$\text{Total Chl Biofilm } (\mu\text{g}) = (\text{OD}_{663} - \text{OD}_{750}) \times 12.7 \times (\text{Dilution Factor}) \times V_{total}$$

- b. To account for the contribution of remaining planktonic cells in the biofilm fraction, subtract the planktonic chlorophyll present in the fraction's volume from the above value:

$$\text{Total Chl Biofilm Corrected } (\mu\text{g}) = (\text{Chl concentration in biofilm}) \times V_{total} - (\text{Chl concentration in planktonic}) \times (\text{Volume of Detached biofilm})$$

3. Total chlorophyll in the culture is the sum of the amount of chlorophyll in planktonic cells and in the biofilm.

Total Chl (μg) = Total Chl Biofilm + Total Chl Planktonic

Or

Total Chl (μg) = Total Chl Biofilm Corrected + Total Chl Planktonic

4. The extent of biofilm formation is presented as percent of chlorophyll in planktonic cells.

Chl in suspension (%) = (Total Chl Planktonic)/(Total Chl) x 100%

Notes

1. All work with cyanobacterial cultures is carried out under sterile conditions using a laminar flow hood.
2. The sponge plug is punched using an awl.
3. For bubbling of a small number of cultures, a mini aquarium pump (e.g., JBL ProAir a50, Art. No. 6054600) may be used instead of a compressor. As an alternative to mixing air and CO₂ one may purchase 3% CO₂ cylinders.
4. A previous study revealed that a biofilm inhibitory substance is present in extracellular fluids from wild type culture as well as in old cultures of the biofilm forming mutant, T2SEQ (Schatz *et al.*, 2013). To avoid possible inhibitor accumulation in cultures of biofilm-forming mutants, cultures are initiated by inoculation of cells grown from solid agar. Make sure to dilute these starter cultures for 'biofilm quantification assays' while yet at exponential phase of growth.
5. Starter cultures are inoculated from cultures grown on solid agar. Due to variations in the inoculum, the resulting liquid cultures vary in cell density. The dilution steps are important to ensure that experiments are initiated at the same density and a similar physiological state.
6. First, take an aliquot from the upper part of the suspended culture and then remove the rest of the suspended fraction to avoid contamination of the extraction aliquot with biofilm cells.

Recipes

1. Liquid BG11-medium (based on Stanier *et al.*, 1971)
 - a. Stock I, 100x concentrated (autoclave)
 - 150.00 g/L NaNO₃
 - 6.50 g/L MgSO₄·7H₂O
 - 3.60 g/L CaCl₂·2H₂O
 - b. Stock II, 100x concentrated (autoclave)
 - 3.05 g/L K₂HPO₄
 - 0.10 g/L Na₂Mg-EDTA

- c. Stock III, 100x concentrated (prepare fresh, no need to sterilize)
 - 0.60 g/L $C_6H_{11}FeNO_7$
 - 0.60 g/L $C_6H_8O_7$
- d. Stock V, 1,000x concentrated (sterilize by filtration using Filter Stericup [0.22 μ m])
 - 2.86 g/L H_3BO_3
 - 1.84 g/L $MnCl_2 \cdot 4H_2O$
 - 0.22 g/L $ZnSO_4 \cdot 7H_2O$
 - 0.39 g/L $NaMoO_4 \cdot 2H_2O$
 - 0.08 g/L $CuSO_4 \cdot 5H_2O$
 - 0.05 g/L $Co(NO_3)_2 \cdot 6H_2O$
- e. HEPES 119.15 g/L (25x concentrated) titrated to pH 8.0 with 10 N NaOH (399.97 g/L) (autoclave)
- f. Preparation of 1 L liquid-BG11:

Start with ~500 ml DDW, add the required volume from the stock solutions and complete with DDW up to 1 L. Solution is sterilized by autoclaving (sterilization time: 40 min). Immediately remove from autoclave once the cycle is over (Long exposure to the high temperature may affect precipitation of minerals). Do not autoclave the medium more than once. Mix before aliquoting the medium to growth tubes

Note: To reproducibly observe biofilm formation, it is crucial to use freshly made stock III each time liquid BG11 is prepared. Autoclaved BG11 is stored at room temperature. The medium should be used within 4 days to inoculate cultures.

2. Solid BG11-medium

- a. Stock solutions as for liquid BG11
- b. Sodium bicarbonate 42 g/L (100x concentrated)
- c. Sodium thiosulfate 470 g/L (100x concentrated)
- d. Prepare two-fold concentrated liquid BG11 (BG11X2) by diluting appropriate volumes from stock solutions
- e. Prepare 3% Bacto agar in DDW
- f. Autoclave BG11X2 and 3% agar in separate containers
- g. Immediately following sterilization, combine equal volumes of BG11X2 and 3% agar to get 1.5% agar in BG11
- h. Cool to ~50-55 °C, dilute 100-fold $NaHCO_3$ and $Na_2S_2O_3$ stock solutions (at this step antibiotics may be added):
 - Spectinomycin dihydrochloride pentahydrate (final concentration: 50 μ g/ml)
 - Kanamycine sulphate monohydrate (final concentration: 50 μ g/ml)
 - Gentamycin sulphate (final concentration: 7.5 μ g/ml)
 - Chloramphenicol (final concentration: 7.5 μ g/ml)
- i. Pour ~50 ml into 90 mm Petri-dish

Acknowledgments

This protocol is modified from previous studies (Schatz *et al.*, 2013; Parnasa *et al.*, 2016). Growth medium is based on the protocol by (Stanier *et al.*, 1971). Chlorophyll measurements are based on previous protocols (Arnon, 1949; Ritchie, 2006). Rakefet Schwarz and Susan Golden are supported by the program of the National Science Foundation and the US-Israel Binational Science Foundation (NSF-BSF 2012823). This study was also supported by a grant from the Israel Science Foundation (ISF 1406/14) to Rakefet Schwarz.

References

1. Agostoni, M., Waters, C. M. and Montgomery, B. L. (2016). [Regulation of biofilm formation and cellular buoyancy through modulating intracellular cyclic di-GMP levels in engineered cyanobacteria](#). *Biotechnol Bioeng* 113(2): 311-319.
2. Arnon, D. I. (1949). [Copper enzymes in isolated chloroplasts. Polyphenoloxidase in beta vulgaris](#). *Plant Physiol* 24(1): 1-15.
3. Enomoto, G., Ni Ni, W., Narikawa, R. and Ikeuchi, M. (2015). [Three cyanobacteriochromes work together to form a light color-sensitive input system for c-di-GMP signaling of cell aggregation](#). *Proc Natl Acad Sci U S A* 112(26): 8082-8087.
4. Enomoto, G., Nomura, R., Shimada, T., Ni Ni, W., Narikawa, R. and Ikeuchi, M. (2014). [Cyanobacteriochrome SesA is a diguanylate cyclase that induces cell aggregation in *Thermosynechococcus*](#). *J Biol Chem* 289(36): 24801-24809.
5. Fisher, M. L., Allen, R., Luo, Y. and Curtiss, R., 3rd (2013). [Export of extracellular polysaccharides modulates adherence of the *Cyanobacterium synechocystis*](#). *PLoS One* 8(9): e74514.
6. Jittawuttipoka, T., Planchon, M., Spalla, O., Benzerara, K., Guyot, F., Cassier-Chauvat, C. and Chauvat, F. (2013). [Multidisciplinary evidences that *Synechocystis* PCC6803 exopolysaccharides operate in cell sedimentation and protection against salt and metal stresses](#). *PLoS One* 8(2): e55564.
7. Nagar, E. and Schwarz, R. (2015). [To be or not to be planktonic? Self-inhibition of biofilm development](#). *Environ Microbiol* 17(5): 1477-1486.
8. Oliveira, P., Pinto, F., Pacheco, C. C., Mota, R. and Tamagnini, P. (2015). [HesF, an exoprotein required for filament adhesion and aggregation in *Anabaena* sp. PCC 7120](#). *Environ Microbiol* 17(5): 1631-1648.
9. O'Toole, G. A. and Kolter, R. (1998). [Initiation of biofilm formation in *Pseudomonas fluorescens* WCS365 proceeds via multiple, convergent signalling pathways: a genetic analysis](#). *Mol Microbiol* 28(3): 449-461.

10. Parnasa, R., Nagar, E., Sendersky, E., Reich, Z., Simkovsky, R., Golden, S. and Schwarz, R. (2016). [Small secreted proteins enable biofilm development in the cyanobacterium *Synechococcus elongatus*](#). *Sci Rep* 6: 32209.
11. Ritchie, R. J. (2006). [Consistent sets of spectrophotometric chlorophyll equations for acetone, methanol and ethanol solvents](#). *Photosynth Res* 89(1): 27-41.
12. Schatz, D., Nagar, E., Sendersky, E., Parnasa, R., Zilberman, S., Carmeli, S., Mastai, Y., Shimoni, E., Klein, E., Yeger, O., Reich, Z. and Schwarz, R. (2013). [Self-suppression of biofilm formation in the cyanobacterium *Synechococcus elongatus*](#). *Environ Microbiol* 15(6): 1786-1794.
13. Schwarzkopf, M., Yoo, Y. C., Huckelhoven, R., Park, Y. M. and Proels, R. K. (2014). [Cyanobacterial phytochrome2 regulates the heterotrophic metabolism and has a function in the heat and high-light stress response](#). *Plant Physiol* 164(4): 2157-2166.
14. Stanier, R. Y., Kunisawa, R., Mandel, M. and Cohen-Bazire, G. (1971). [Purification and properties of unicellular blue-green algae \(order *Chroococcales*\)](#). *Bacteriol Rev* 35(2): 171-205.
15. Zavřel, T., Sinetova, M.A., and Červený, J. (2015). [Measurement of chlorophyll a and carotenoids concentration in *Cyanobacteria*](#). *Bio-protocol* 5: 1-5.